



Veraview iX

Intraoral Imaging

- Accurate LED Positioning Guide
- High Frequency DC Generator Reduces Patient Dosage
- Optimized Time Settings for Film and Digital Exposures
- Simplified Operation with Illuminated Pictogram Menu
- Available in two versions: wall-mounted or with patient's chair

Thinking ahead. Focused on life.


Quick and Accurate Alignment by LED Positioning Guide



iX-R model with patient's chair

■ LED Positioning Guide

- Two LED light beams converge into a single point for visual confirmation of correct alignment for accurate X-rays.
- The cone is short and easy to manipulate which also facilitates correct positioning. It is also transparent for superior visibility.
- The lightweight head and arm can easily be maneuvered with a delicate touch.

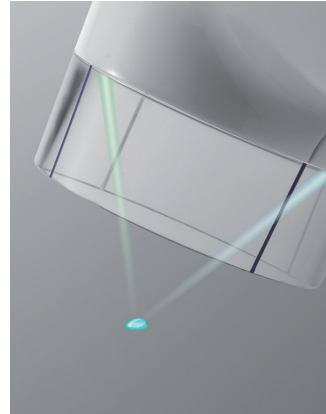

■ Constant Voltage DC Generator

- By utilizing a high frequency, direct current (DC) generator, this unit produces 'hard' X-rays which offer reduced dosage compared to units that operate with an alternating current.
- Tube voltage can be set at 60 kV or 70 kV.

■ Easy Seat Height Adjustment (Type R)


Seat height is easily adjusted by pulling up the lever for the gas cylinder under the seat.

■ Illuminated Pictogram Menu

Menu items illuminate to confirm selections and help prevent mistakes.

■ Exposure Settings

Digital or film scan options are available. An optimized exposure time is set after parameter selections are made.

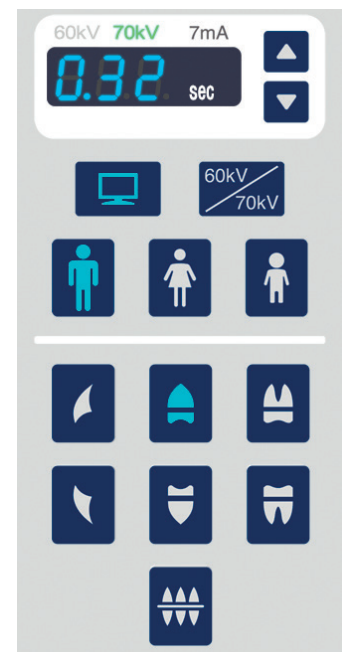

Order no.

1240 Verview iX Type R



iX-WA wall-mounted version

Order no.

1220 Verview iX Type WA550

1230 Verview iX Type WA800

1235 Verview iX Type WA1000

Specifications

Model	Verview iX (wall-mounted version)	Verview iX-R (version with patient's chair)
Voltage	100-120 / 220-240 VAC	100-120 / 220-240 VAC
Power Consumption	1,9 kVA	1,9 kVA
Tube Potential	60 / 70 kV	60 / 70 kV
Tube Current 7 mA	7 mA	7 mA
Exposure Time	0.01~2 seconds	0.01~2 seconds
Weight	approx. 23 kg	approx. 25 kg

Development and Manufacturing

J. Morita Mfg. Corp.
680 Higashihama Minami-cho, Fushimi-ku
Kyoto 612-8533, Japan
T +81. 75. 611 2141, F +81. 75. 622 4595
www.morita.com

Distribution

J. Morita Corporation
33-18, 3-Chome, Tarumi-cho Suita City, Osaka, 564-8650 Japan
T +81. 6. 6380 1521, F +81. 6. 6380 0585
www.morita.com/asia

Subject to technical changes and errors.
JME JB16-005EN 0216 *1

J. Morita Europe GmbH

Justus-von-Liebig-Str. 27a, 63128 Dietzenbach, Germany
T +49. 6074. 836 0, F +49. 6074. 836 299
www.morita.com/europe

J. Morita Corporation Australia & New Zealand

247 Coward Street, Mascot, NSW 2020, Australia
T +61. 2. 9667 3555, F +61. 2. 9667 3577
www.morita.com/anz

J. Morita Middle East

Saraya Al Wessam Building, Appt 902 Cross
Roads Tag Al Roasaa & Abu Heif SabaPacha, Alexandria, 21311, Egypt
T +20. 3. 58 222 94, F +20. 3. 58 222 96
www.morita.com/middleeast